


ToeGrip[®] Classic

First Hammer Toe
PEEK implanted
with over 50'000
units implanted since
2010

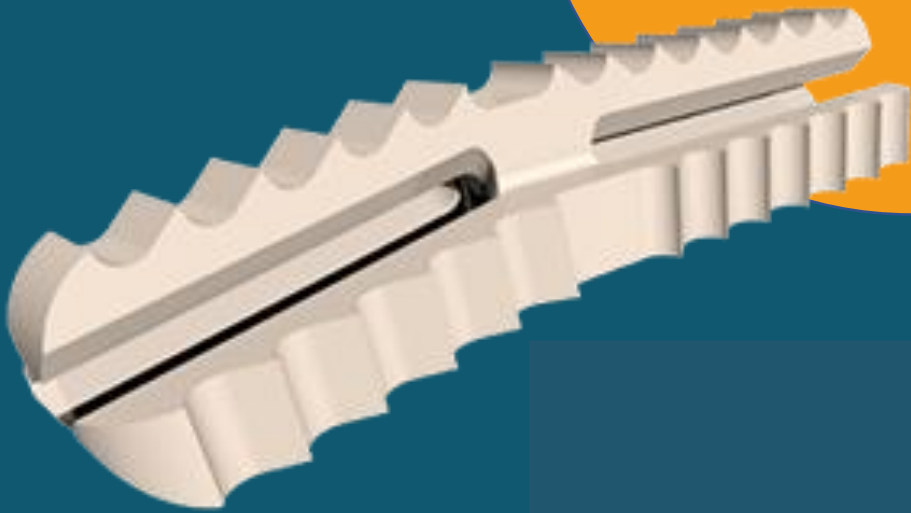

SynchroMedical

The interphalangeal reference

ToeGrip® Classic


Interphalangeal Implant

DESCRIPTION

The ToeGrip® Classic implant consists of a single part implant with three flexible intramedullary prongs inserted into each part of the phalanx. The prongs are attached by a T-shaped structure. The attachment concept is based on a press-fit contact due to its tapered shape and the macrostructures solidly anchored once impacted, thus avoiding any displacement of the implant.

The ToeGrip® Classic implant is intended for single use only and is available in a range of 5 sizes with 3 possible angulations: 0°, 10° or 20°.

ToeGrip® Classic is made of PEEK as per standard ASTM F2026.

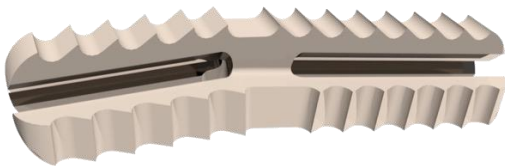

Clinical Evidence

The ToeGrip® Classic has over 10 years of clinical experience with publications in peer reviewed journals. Since 2010 over 50'000 cases have been successfully treated.

Publications

Interphalangeal Arthrodesis of the Toe With a New Radiolucent Intramedullary Implant.

Foot & Ankle Specialist Review. Dec 2015

Diebold P, Rocher H, Determe P, Cermolacce C, Guillo S, Averous C, Lieber Wackenhein F.

A Midterm Review of Lesser Toe Arthrodesis With an Intramedullary Implant.

Foot & Ankle Specialist Review. Oct 2017

James Lee Harmer, FCPodS, MSc, BSc(Hons), Anthony Wilkinson, FCPodS, MSc, DpodM, and Anthony John Maher, FCPodS, MSc, BSc(Hons).

ORDERING INFORMATION

Instruments

TG0240-Starter		P1 COMPACTOR - STARTER - AO COUPLING
TG0240-XS		P1 COMPACTOR - XS - AO COUPLING
TG0240-SS		P1 COMPACTOR - S - AO COUPLING
TG0240-MM		P1 COMPACTOR - M - AO COUPLING
TG0240-LL		P1 COMPACTOR - L - AO COUPLING

TG0241-Starter		P2 COMPACTOR - STARTER - AO COUPLING
TG0241-XS		P2 COMPACTOR - XS - AO COUPLING
TG0241-SS		P2 COMPACTOR - S - AO COUPLING
TG0241-MM		P2 COMPACTOR - M - AO COUPLING
TG0241-LL		P2 COMPACTOR - L - AO COUPLING

TG0203		REDUCER
TG0009		HOLDER SIZE XS
TG0006		HOLDER SIZE S
TG0007		HOLDER SIZE M
TG0008		HOLDER SIZE L

TG0010		CALCAR REAMER
--------	--	---------------

HV9010		AO HANDLE
--------	--	-----------

TG1013		ToeGrip® INSTRUMENT TRAY
--------	--	--------------------------

Implants

TGXS00	ToeGrip® Classic EXTRA-SMALL 0°
TGXS10	ToeGrip® Classic EXTRA-SMALL 10°
TGXS20	ToeGrip® Classic EXTRA-SMALL 20°
TGSS00	ToeGrip® Classic SMALL 0°
TGSS10	ToeGrip® Classic SMALL 10°
TGSS20	ToeGrip® Classic SMALL 20°
TGMM00	ToeGrip® Classic MEDIUM 0°
TGMM10	ToeGrip® Classic MEDIUM 10°
TGMM20	ToeGrip® Classic MEDIUM 20°
TGLL00	ToeGrip® Classic LARGE 0°
TGLL10	ToeGrip® Classic LARGE 10°
TGLL20	ToeGrip® Classic LARGE 20°