


ToeGrip® Evo

Anatomically
adapted based on
the experience of
over 50'000
ToeGrip® Classic

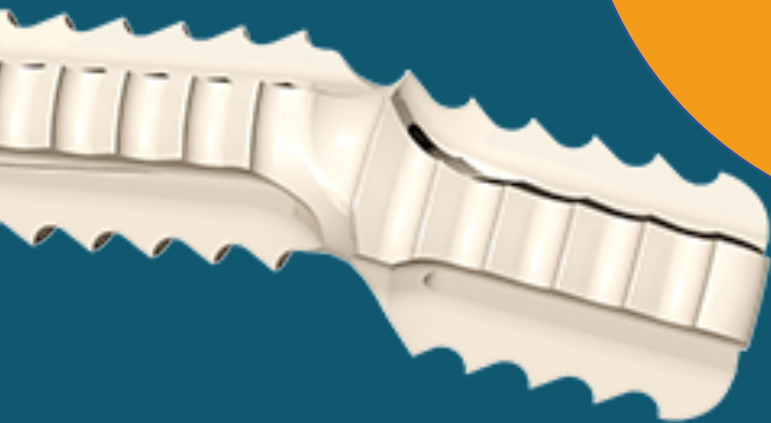

SynchroMedical

The interphalangeal reference

ToeGrip® Evo

Interphalangeal implant

DESCRIPTION


The ToeGrip® EVO implant consists of a single part implant with three flexible intramedullary prongs on each side. The prongs are inserted into each the proximal and distal phalanx.

The prongs are connected by a cylindrical structure with a sagittal height offset between the proximal and distal phalanx to mimic the toe anatomy.

The device anchorage relies on a press-fit contact that is supported by its tapered shape and the macrostructure, which solidly anchored once impacted.

Primary stability and avoiding any displacement are key points for a successful bone consolidation.

ToeGrip® EVO is made from PEEK as per standard ASTM F2026.

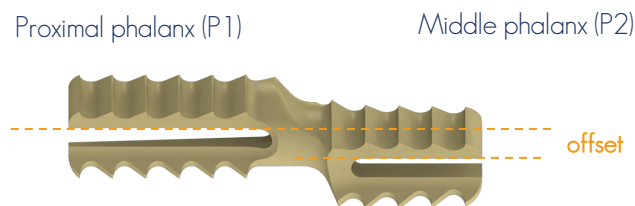

Clinical Evidence

The ToeGrip® Classic has over 10 years of clinical experience with publications in peer reviewed journals. Since 2010 over 50'000 cases have been successfully treated.

Publications

Interphalangeal Arthrodesis of the Toe With a New Radiolucent Intramedullary Implant.

Foot & Ankle Specialist Review. Dec 2015

Diebold P, Rocher H, Determe P, Cermolacce C, Guillo S, Averous C, Lieber Wackenheim F.

A Midterm Review of Lesser Toe Arthrodesis With an Intramedullary Implant.

Foot & Ankle Specialist Review. Oct 2017

James Lee Harmer, FCPodS, MSc, BSc(Hons), Anthony Wilkinson, FCPodS, MSc, DpodM, and Anthony John Maher, FCPodS, MSc, BSc(Hons).

ORDERING INFORMATION

Implants

TGEXS00	ToeGrip® EVO EXTRA-SMALL 0°
TGEXS10	ToeGrip® EVO EXTRA-SMALL 10°
TGESS00	ToeGrip® EVO SMALL 0°
TGESS10	ToeGrip® EVO SMALL 10°
TGEMM00	ToeGrip® EVO MEDIUM 0°
TGEMM10	ToeGrip® EVO MEDIUM 10°

Instruments

	TG0240-Starter	P1 COMPACTOR - STARTER
	TG0240-XS	P1 COMPACTOR - XS
	TG0240-SS	P1 COMPACTOR - S
	TG0240-MM	P1 COMPACTOR - M
	TG0010	CALCAR REAMER
	HV9010	AO HANDLE
	TG0203	REDUCER
	TG0020	HOLDER SIZE XS 0°
	TG0025	HOLDER SIZE XS 10°
	TG0021	HOLDER SIZE S 0°
	TG0026	HOLDER SIZE S 10°
	TG0022	HOLDER SIZE M 0°
	TG0027	HOLDER SIZE M 10°
	TGE9010	ToeGrip® EVO INSTRUMENT TRAY